

KEN WASHER DISINFECTOR IQ6

Intelligent Quality

- The Future of Efficient Disinfection

- ✓ Easy to use
- ✓ Easy to service
- ✓ Easy to choose

The IQ series is a modular built range that meets all relevant standards such as ISO 15883, UL to name a few.

The elegant design provides a smooth surface in tempered glass that is easy to clean. The IQ6 has downward sliding doors that allows for manual loading, or the use of a unique automatic loading that provides high flexibility and uses very little space.

The IQ6 is easy to operate by touch panel with swipe function, showing all relevant information of the program needed by the operators.


IQ6 has the lowest energy and water consumption in the market, to keep operating costs at a minimum and contribute to a clean environment.

IQ6 comes with a wide range of racks and modules and holds up to 18 DIN baskets.

All data from the IQ6 are fully exchangeable with common computer systems.



HEALTHCARE

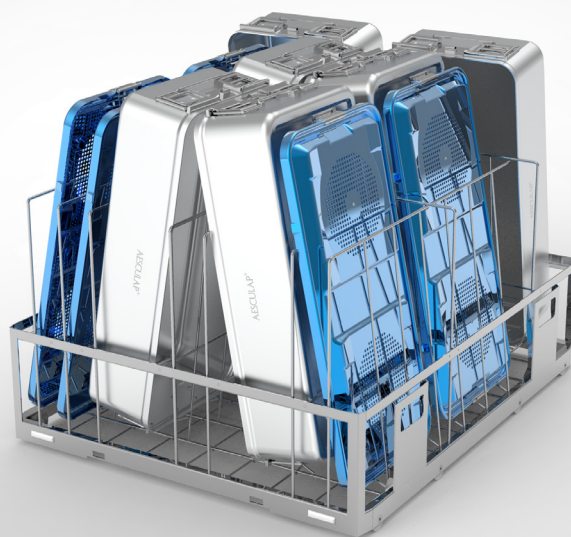

A wide range of racks, modules and even special demands equipment can be provided to solve your disinfection needs.

Advantages:

Environment

- Low energy costs
- Low noise level
- Low chemical usage
- Very low water consumption

Technology

- RFID recognition
- Robotic loading / unloading
- Service access from the front

Hygiene safety

- EN 62366 tested
- EN ISO 15883-1-2-5

Validation of all processes

- Open programs structure
- Full process validation
- Exchange of data

Description	Values
Outer dimensions	H: 1985 x W: 820 x D: 936 mm
Operative chamber dimensions	H: 690 x W: 705 x D: 850 mm
Weight excl./incl. packaging	340 / 405 kg
Operative volume, wash chamber	413 L
Total volume, wash chamber	492 L
Capacity DIN-baskets 480x250x50	18
Cycle time, standard	35-45 min (fast cycle machine < 30 min)
Exhaust connections	ø100 mm 385 m3 per hour
Water connections	2 x 3/4 (3rd is optional)
Cold water	200-800 kPa, ≥ 15 L/min, 5-30°C, 0-30° dH
Hot water	200-800 kPa, ≥ 15 L/min, 30-70°C, 0-3° dH
Demineralized / RO water	200-800 kPa, ≥ 15 L/min, 5-70°C - option
Washing / Drying	<56.5 dB (A) / <60.5 dB (A)
Water consumption per phase	Approx. 22-25 L (standard 5-level rack)
Power supply options	3x200-480V ± N + PE, 50/60 Hz
Connected load	14-26 kW, 40A
Programs available	10 standard + 30 optional
Drain	ø50 mm

Contact a KEN partner

Clean, Lean & Green